


Promotor[®]

Nucleic Acid Extraction Solution

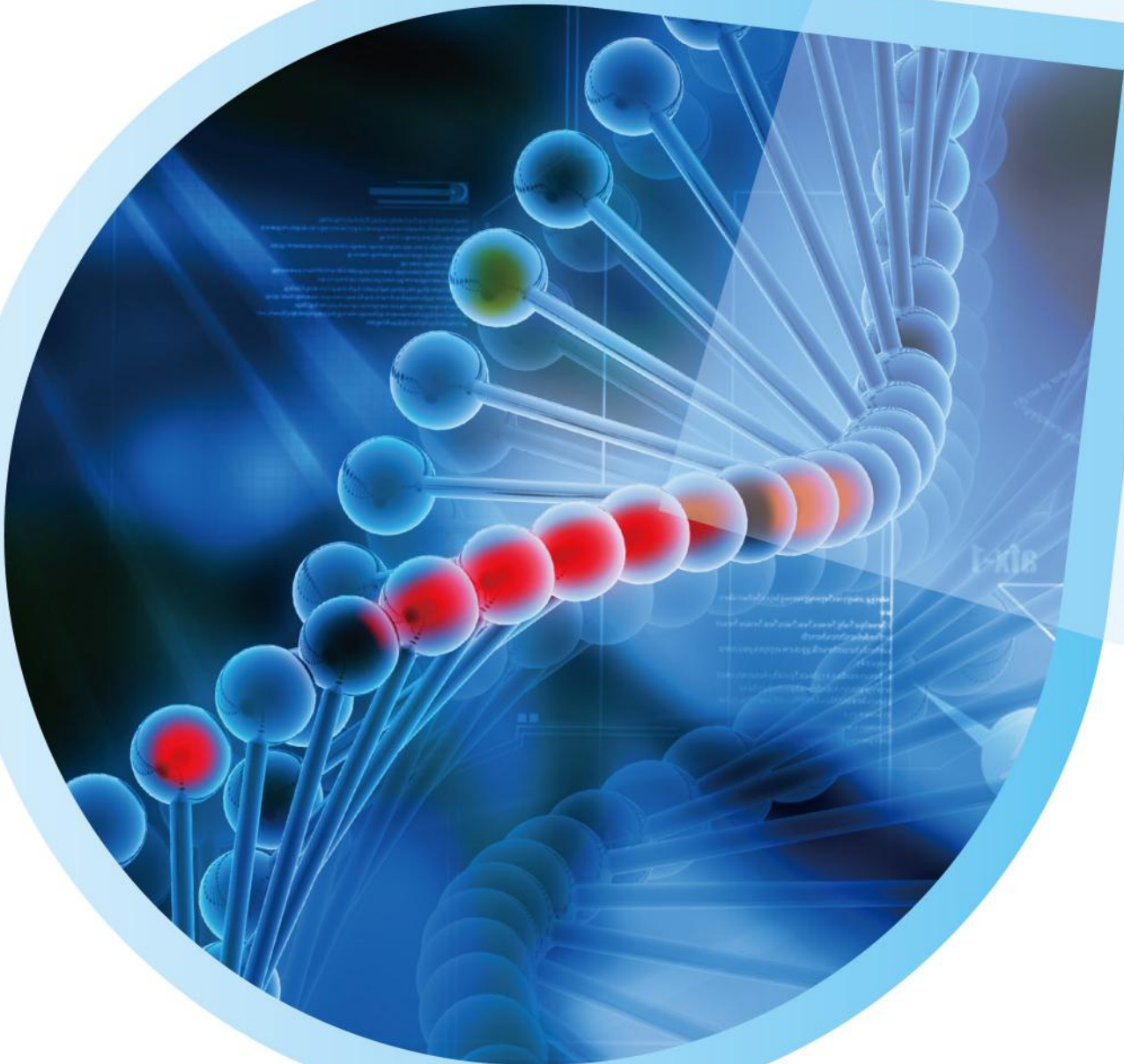

ACON[®]

►►► NES-32 Nucleic Acid Extraction System



Features & Benefits

- Up to 32 samples/run
- High purity & high yield
- Ready-to-use reagents
- 8-inch large color touch screen
- Aviation material heating base
- UV lamp and disposable consumables eliminate contamination



DNA


RNA

Genomic DNA

DNA/RNA



Principle

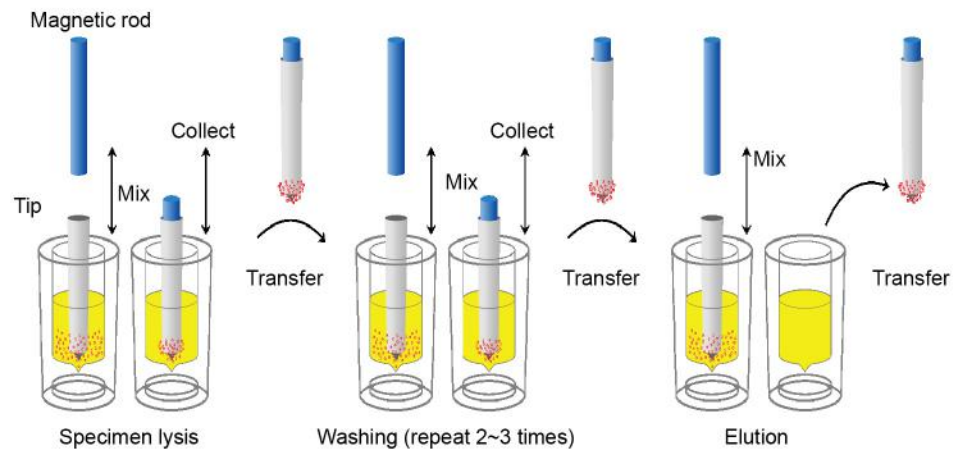

Specifications

Parameters	Specifications
Method	Magnetic Beads
Processing Volume	20 to 1,000µL
Sample Capacity	Up to 32 samples per run
Sample Type	Serum, Plasma, Swab, Feces, Urine
Well-to-Well CV	≤1.5%
Collection Rate of Magnetic Particles	≥98%
Heating Temperature	Lysis Temp (from RT to 120°C) Elution Temp (from RT to 120 °C)
Sterilization	UV Lamp
Processing Time	12-40 min
Dimensions	437×512×470 mm (W×D×H)
Instrument Port	USB/LAN
Power Source	100-240 VAC, 50 Hz/60 Hz, 600VA

▶▶▶ Viral Nucleic Acid Isolation Kits(Spin Column)



RNA



DNA/RNA

Use a simple, rapid protocol to generate high-quality templates for a broad range of downstream applications.

- Fast procedures and easy to use.
- Isolation and purification of high-quality viral DNA/RNA.
- High efficient recovery yield.
- Pure DNA/RNA ready for downstream applications.
- No phenol-chloroform extraction.



High Pure Viral Nucleic Acid Isolation Kits Workflow

1 Lysis



Specimen lysis

2 Bind



Load Specimen

3 Wash



Washing and dry

4 Elute



Nucleic Acid elution

Specifications

Parameters	Specifications
Format	Spin Column
Technology	Silica-membrane Technology
Sample Volume	200 μ L
Elution Volume	50 μ L
Sample Type	Serum, Plasma, Swab
Yield	>90% recovery
Processing	Manual (Centrifugation)
Processing Time per run	15-30 min
Purification of total RNA, miRNA, poly A+ Mrna, DNA or protein	Viral DNA or Viral RNA

Ordering Information

Product Name	Catalog No.	Components		Type	Specimen	Kit Size
NES-32 Nucleic Acid Extraction System	P111-1011√*	Analyzer Manual		Magnetic Separation Analyzer	Serum, Plasma, Whole blood, Swab, Feces or Urine	1 Unit/Kit
Nucleic Acid Extraction Kit	P121-1301√*	Extraction Plate Protease K	Carrier RNA Tip Comb	Magnetic Beads	Serum, Plasma, Swab, Sputum	32 Tests/Kit
Nucleic Acid (DNA) Extraction Kit	P121-1101√*	Tip Combs 96 Deep Well Plates Magnetic Beads Liquid Paraffin Proteinase K	Wash Solution-1 Wash Solution-2 Lysis Solution Elution Solution	Magnetic Beads	Serum, Plasma, Swab, Feces or Urine	32 Tests/Kit
Whole Blood Genomic DNA Extraction Kit	P121-1111√*	Tip Combs 96 Deep Well Plates Magnetic Beads Liquid Paraffin Proteinase K	Wash Solution-1 Wash Solution-2 Lysis Solution Elution Solution	Magnetic Beads	Whole blood	32 Tests/Kit
Nucleic Acid (RNA) Extraction Kit	P121-1121√*	Tip Combs 96 Deep Well Plates Magnetic Beads Liquid Paraffin	Wash Solution-1 Wash Solution-2 Lysis Solution Elution Solution	Magnetic Beads	Serum, Plasma, Swab, Feces or Urine	32 Tests/Kit
Nucleic Acid (DNA) Extraction Kit (ready-to-use)	P121-1201√*	Tip Combs Extraction Plate Proteinase K		Magnetic Beads	Serum, Plasma, Swab, Feces or Urine	32 Tests/Kit
Whole Blood Genomic DNA Extraction Kit (ready-to-use)	P121-1211√*	Tip Combs Extraction Plate Proteinase K		Magnetic Beads	Whole blood	32 Tests/Kit
Nucleic Acid (RNA) Extraction Kit (ready-to-use)	P121-1221√*	Tip Combs Extraction Plate		Magnetic Beads	Serum, Plasma, Swab, Feces or Urine	32 Tests/Kit
Viral DNA/RNA Isolation Kit	P121-1291√*	Lysis Solution Wash Solution-1 Wash Solution-2 Elution Solution	Proteinase K Spin Columns Collection Tubes	Spin Column	Serum, Plasma or Swab	32 Tests/Kit 50 Tests/Kit
Viral Nucleic Acid Isolation Kit	P121-1051√*	Kit 1: Buffer L Buffer WA Buffer WB	Buffer TE Collection Tubes (2mL) Spin Columns	Spin Column	Serum, Plasma or Swab	32 Tests/Kit
		Kit 2: Carrier RNA				

√ CE marked *FDA Device Listing

We also offer other diagnostic and medical products for:
Blood Glucose, Lateral Flow, Urinalysis Products, POCT Products, Immunoassay Products and more.
 Contact us for worldwide distribution and custom manufacturing (OEM) opportunities.

✓ CE Marked for sale in the European Community 

All right reserved. Programs and Specifications subject to changes without prior notice, please contact local distributors.



aconlabs.com

ACON Laboratories, Inc.
 5850 Oberlin Drive, #340, San Diego, CA 92121, USA
 Tel: 1-858-875-8000
 Fax: 1-858-200-0729
 E-mail : info@aconlabs.com